


YORLAB[®] SPINMAT CYTO CONSUMABLES **Model No.: YSL / YS-CC101**



**CYTO FUNNELS WITH
FILTER PAPERS**



**CYTO FUNNELS
ONLY**

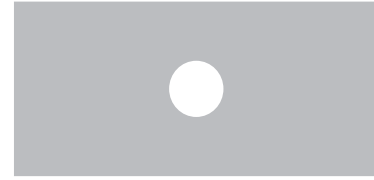

**FILTER PAPER
ONLY**

Specifications:

- Development for single-use with the Shandon & other Similar Cyto centrifuges and other cytology slide preparation systems of centrifugation mechanism on the market, the YORLAB SpinMAT cyto consumables are much safer than re-usable cyto chambers and lower the risk of contamination.
- Compared with most of the other similar products on the market, the YORLAB SpinMAT Cyto consumable are more trustable depending our serious and professional R & D and production.

The YORK SpinMAT Cyto consumables are supplied in :

- Filter Paper Only .
- Cyto Funnels Only.
- Cyto Funnels w /Filter papers.
- SS Slide Clip Holding.
- YORK SpinMat Cyto Funnel only 50/PK
- YORK SpinMat Cyto Funnel W/White Filter Paper 50/PK
- YORK SpinMat Cyto Funnel w/Tan Filter paper 50/PK
- YORK SpinMat White Filter Paper only 200/PK
- YORK SpinMat Tan Filter Paper only 200/PK
- Stainless Steel Clip Holding Disposable
- Funnel Sample Chamber
- SS Slide clip for Holding cytology disposable funnel sample chamber against the microscope slide for preparation in the Cyto centrifuge.
- Eliminates contact between contaminated SS Slide clip

STAINLESS STEEL CLIP HOLDING



for holding disposable funnels and laboratory personnel.

- Clips are autoclavable and reusable.

Description :

- Stainless Steel Slide Clip.

How the Single Chamber Funnel Works

- The Single Chamber Funnel fits snugly against a specially marked slide and is held in place by a Stainless Steel Slide Clip. Twelve Stainless Steel Slide Clip assemblies fit into the Cyto centrifuge head.
- After cyto centrifugation, the Stainless Steel Slide clip easily unlocks to dispose of the sample chamber and to remove the slide.